

Date : 2021-09-30
No. : HC21060354

Page 1 of 8

Applicant(Code:01325900) : Savewo Limited
1/F
266-270 Texaco Road
Tsuen Wan NT HK

Description of Sample(s) : One submitted sample said to be Savewo 3DMask Ultra.
Country of Origin : Hong Kong

Sample(s) Received Condition(s): In intact original package under ambient temperature

Date Sample(s) Received : 2021-06-08

Date Tested : 2021-09-15 to 2021-09-29

Investigation Requested : Performance Test as per ASTM F2100-19
1. Bacterial Filtration Efficiency (BFE) %
– *Staphylococcus aureus* (ATCC 6538)
2. Particulate Filtration Efficiency (PFE) %
3. Differential Pressure
4. Synthetic Blood Penetration
5. Flammability to Class 1

Remark(s) : Tests were conducted on the sample after accelerated aging at 60°C for 97 days according to ASTM F1980.



LAU Yuk Kuen, Joey
Authorized Signatory





STC
Test Report

Date : 2021-09-30
No. : HC21060354

Page 2 of 8

Requirement:

| Performance Test as per ASTM F2100-19 | Level 1 | Level 2 | Level 3 |
|---|--|---|----------|
| Bacterial Filtration Efficiency (BFE) % – <i>Staphylococcus aureus</i> (ATCC 6538) | ≥95% | ≥98% | |
| Particulate Filtration Efficiency (PFE) % | ≥95% | ≥98% | |
| Differential Pressure (ΔP) | <5.0 mmH ₂ O/cm ² | <6.0 mmH ₂ O/cm ² | |
| Resistance to Penetration by Synthetic Blood | 80 mmHg | 120 mmHg | 160 mmHg |
| Flame Speed (Flammability to Class 1) | Class 1 (The time of flame spread is 3.5 seconds or more) | | |

Summary:

| Performance Test as per ASTM F2100-19 | Savewo 3DMask Ultra |
|--|---------------------|
| | Level 3 |
| Bacterial Filtration Efficiency(BFE) % – <i>Staphylococcus aureus</i> (ATCC 6538) | Pass |
| Particulate Filtration Efficiency (PFE) % | Pass |
| Differential Pressure (ΔP) | Pass |
| Resistance to Penetration by Synthetic Blood Penetration | Pass |
| Flame Speed (Flammability to Class 1) | Pass |

Note: An acceptable quality limit of 4% shall be used for all required testing to establish conformance of medical face masks to a specific performance class.

Date : 2021-09-30
No. : HC21060354

Page 3 of 8

Test Result(s):

1. Bacterial Filtration Efficiency (BFE) %

Test method: ASTM F2100-19 9.1 & ASTM F2101-19

Summary: The BFE test is performed to determine the filtration efficiency of test articles by comparing the bacterial control counts upstream of the test article to the bacterial counts downstream. A suspension of *Staphylococcus aureus* was aerosolized using a nebulizer and delivered to the test article at a constant flow rate and fixed air pressure. The challenge delivery was maintained at $1.7 - 3.0 \times 10^3$ colony forming units (CFU) with a mean particle size (MPS) of $3.0 \pm 0.3\mu\text{m}$. The aerosols were drawn through a six-stage, viable particle, Andersen sampler for collection. This test method complies with ASTM F2101-19.

All test method acceptance criteria were met.

| Specimen(s) | Savewo 3DMask Ultra |
|-------------|---------------------|
| 1 | 99.6% |
| 2 | 99.6% |
| 3 | 99.7% |
| 4 | 99.8% |
| 5 | 99.9% |

- Notes :
- Challenge bacteria : *Staphylococcus aureus* (ATCC 6538)
 - Positive control average : 1724 CFU
 - Negative control average : <1 CFU
 - Mean particle size : 3.1 μm
 - Testing side : Outside of specimen
 - Testing area : 49 cm²
 - Precondition : Minimum of 4 hours at (21 \pm 5) °C and (85 \pm 5) % relative humidity (RH)

Date : 2021-09-30
No. : HC21060354

Page 4 of 8

2. Particulate Filtration Efficiency (PFE) %

Test method: ASTM F2100-19 9.3 & ASTM F2299-17

Summary: This procedure was performed to evaluate the non-viable particle filtration efficiency (PFE) of the test article. Monodispersed polystyrene latex spheres (PSL) were nebulized (atomized), dried, and passed through the test article. The particles that passed through the test article were enumerated using a laser particle counter.

The upstream and downstream particle counts at each position were sampled and recorded. The filtration efficiency was calculated using the average number of particles penetrating the test article (downstream particle count) compared to the average of the upstream particle count.

The procedure employed the basic particle filtration method described in ASTM F2299-17. All test method acceptance criteria were met.

| Specimen(s) | Savewo 3DMask Ultra | | | |
|-------------|-------------------------|---------------------------|---------------------------------|-------|
| | Upstream particle count | Downstream particle count | Resistances to Ventilation (Pa) | PFE % |
| 1 | 99490 | 350 | 90 | 99.6 |
| 2 | 100890 | 360 | 95 | 99.6 |
| 3 | 100580 | 1240 | 90 | 98.8 |
| 4 | 94290 | 930 | 91 | 99.0 |
| 5 | 98580 | 1400 | 88 | 98.6 |

Notes : - Flow rate : 28.3 Litre/min
 - Challenge particles : 0.1 µm PSL
 - Testing area : 50 cm²
 - Testing side : Outside of specimen
 - Testing condition : 18 - 24 °C, 25 -55 % Relative humidity

Date : 2021-09-30
No. : HC21060354

Page 5 of 8

3. Differential Pressure

Test method: ASTM F2100-19 9.2 & EN 14683:2019 + AC:2019, Annex C

Summary: The Differential Pressure test is performed to determine the breathability of test articles by measuring the differential air pressure on either side of the test article using a manometer, at a constant flow rate. This test complies with EN14683:2019 + AC:2019, Annex C and ASTM F2100-19.

All test method acceptance criteria were met.

Sample : Savewo 3DMask Ultra

| Specimen(s) | Test area (in Pa/cm ²) | | | | | Average | |
|-------------|------------------------------------|------|------|------|------|--------------------|------------------------------------|
| | 1 | 2 | 3 | 4 | 5 | Pa/cm ² | mmH ₂ O/cm ² |
| 1 | 45.5 | 48.0 | 48.8 | 41.7 | 49.4 | 46.7 | 4.8 |
| 2 | 49.1 | 43.9 | 48.8 | 43.8 | 46.1 | 46.3 | 4.7 |
| 3 | 47.7 | 48.0 | 49.2 | 48.9 | 47.7 | 48.3 | 4.9 |
| 4 | 47.2 | 45.3 | 55.1 | 50.0 | 50.5 | 49.6 | 5.1 |
| 5 | 45.1 | 46.1 | 45.4 | 49.8 | 50.2 | 47.3 | 4.8 |

Notes : - 1 mmH₂O/cm² = 9.8 Pa/cm²
 - Flow rate: 8 Litre/min
 - Precondition : Minimum of 4 hours at (21±5) °C and (85±5) % relative humidity (RH)

Date : 2021-09-30
No. : HC21060354

Page 6 of 8

4. Synthetic Blood Penetration

Test method: ASTM F2100-19 9.4 & ASTM F1862-17

Summary: This procedure was performed to evaluate surgical facemasks and other types of protective clothing materials designed to protect against fluid penetration. The purpose of this procedure is to simulate an arterial spray and evaluate the effectiveness of the test article in protecting the user from possible exposure to blood and other body fluids. The distance from the target area surface to the tip of the cannula is 30.5cm. A test volume of 2 mL of synthetic blood was employed using the targeting plate method.

This test method was designed to comply with ASTM F1862-17.

Test Pressure: 160mmHg

| Specimen Number | Sawewo 3DMask Ultra |
|------------------------|---------------------|
| 1-20, 22, 24-26, 28-32 | None Seen |
| 21, 23, 27 | Seen |

Requirement:

An acceptable quality limit of 4.0% is met for a normal single sampling plan when ≥ 29 of 32 test specimens show passing result (none seen)

Notes : - Test Side: Outside
 - Precondition : Minimum of 4 hours at (21 ± 5) °C and (85 ± 5) % relative humidity (RH)
 - Testing condition: 18 - 24 °C, 25 -55 % Relative humidity

Date : 2021-09-30
No. : HC21060354

Page 7 of 8

5. Flammability

Test method: ASTM F2100-19 9.5, 16 CFR 1610

Summary: This procedure was performed to evaluate the flammability of plain surface clothing textiles by measuring the ease of ignition and the speed of flame spread. The parameter of time is used to separate materials into different classes, thereby assisting in a judgment of fabric suitability for clothing and protective clothing material. The test procedure was performed in accordance with the test method outlined in 16 CFR Part 1610(a) Step 1 – testing in the original state. Step 2 – Refurbishing and testing after refurbishing, was not performed. All test method acceptance criteria were met.

| Specimen(s) | Savewo 3DMask Ultra | Class |
|-------------|--|-------|
| | Time of spread of flame (Original state) | |
| 1 | Did not ignite | 1 |
| 2 | Did not ignite | |
| 3 | Did not ignite | |
| 4 | Did not ignite | |
| 5 | Did not ignite | |

Notes : - Test Side: Outside
- Orientation: Cross


STC

Test Report

Date : 2021-09-30
No. : HC21060354

Page 8 of 8

Photo(s):

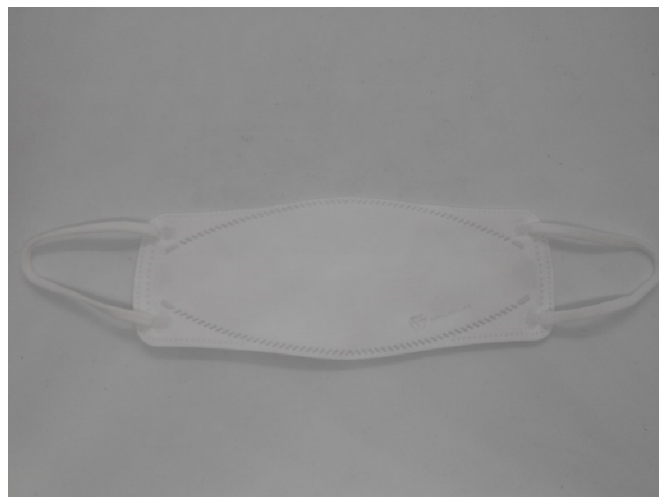

***** End of Test Report *****

Conditions of Issuance of Test Reports

1. All samples and goods are accepted by The Hong Kong Standards & Testing Centre Limited (the "Company") solely for testing and reporting in accordance with the following terms and conditions. The Company provides its services on the basis that such terms and conditions constitute express agreement between the Company and any person, firm or company requesting its services (the "Clients").
2. Any report issued by the Company as a result of this application for testing service (the "Report") shall be issued in confidence to the Clients and the Report will be strictly treated as such by the Company. It may not be reproduced either in its entirety or in part and it may not be used for advertising or other unauthorized purposes without the written consent of the Company. The Clients to whom the Report is issued may, however, show or send it, or a certified copy thereof prepared by the Company to his customer, supplier or other persons directly concerned. Subject to clause 3, the Company will not, without the consent of the Clients, enter into any discussion or correspondence with any third party concerning the contents of the Report, unless required by the relevant governmental authorities, laws or court orders.
3. The Company shall be at liberty to disclose the testing-related documents and/or files anytime to any third-party accreditation and/or recognition bodies for audit or other related purposes. No liabilities whatsoever shall attach to the Company's act of disclosure.
4. The Company shall not be called or be liable to be called to give evidence or testimony on the Report in a court of law without its prior written consent, unless required by the relevant governmental authorities, laws or court orders.
5. The results in Report apply only to the sample as received and do not apply to the bulk, unless the sampling has been carried out by the Company and is stated as such in the Report.
6. When a statement of conformity to a specification or standard is provided, the ILAC-G8 Guidance document (and/or IEC Guide 115 in the electrotechnical sector) will be adopted as a decision rule for the determination of conformity unless it is inherent in the requested specification or standard, or otherwise specified in the Report.
7. In the event of the improper use the report as determined by the Company, the Company reserves the right to withdraw it, and to adopt any other additional remedies which may be appropriate.
8. Sample submitted for testing are accepted on the understanding that the Report issued cannot form the basis of, or be the instrument for, any legal action against the Company.
9. The Company will not be liable for or accept responsibility for any loss or damage howsoever arising from the use of information contained in any of its Reports or in any communication whatsoever about its said tests or investigations.
10. Clients wishing to use the Report in court proceedings or arbitration shall inform the Company to that effect prior to submitting the sample for testing.
11. Subject to the variable length of retention time for test data and report stored hereinto as to otherwise specifically required by individual accreditation authorities, the Company will only keep the supporting test data and information of this test report for a period of three years. The data and information will be disposed of after the aforementioned retention period has elapsed. Under no circumstances shall we provide any data and information which has been disposed of after the retention period. Under no circumstances shall we be liable for damages of any kind, including (but not limited to) compensatory damages, lost profits, lost data, or any form of special, incidental, indirect, consequential or punitive damages of any kind, whether based on breach of contract of warranty, tort (including negligence), product liability or otherwise, even if we are informed in advance of the possibility of such damages.
12. Issuance records of the Report are available on the internet at www.stc.group. Further enquiry of validity or verification of the Reports should be addressed to the Company.

Date : 2021-09-30
No. : HC21060355

Page 1 of 7

Applicant(Code:01325900) : Savewo Limited
1/F
266-270 Texaco Road
Tsuen Wan NT HK

Description of Sample(s) : One submitted sample said to be Savewo 3DMask Ultra.
Country of Origin : Hong Kong

Sample(s) Received Condition(s): In intact original package under ambient temperature

Date Sample(s) Received : 2021-06-08

Date Tested : 2021-09-15 to 2021-09-29

Investigation Requested : Performance Test as per EN 14683:2019 + AC:2019
1. Bacterial Filtration Efficiency (BFE) %
– *Staphylococcus aureus* (ATCC 6538)
2. Differential pressure
3. Microbial cleanliness (Bioburden)
4. Splash resistance

Remark(s) : Tests were conducted on the sample after accelerated aging at 60°C for 97 days according to ASTM F1980.



LAU Yuk Kuen, Joey
Authorized Signatory





STC
Test Report

Date : 2021-09-30
No. : HC21060355

Page 2 of 7

Requirement

| Performance Test as per EN 14683:2019 + AC:2019 | Type I | Type II | Type IIR |
|--|------------------------|---------|------------------------------------|
| Bacterial Filtration Efficiency(BFE) % – <i>Staphylococcus aureus</i> (ATCC 6538) | ≥95% | ≥98% | |
| Breathability (Differential pressure) | <40 Pa/cm ² | | <60 Pa/cm ² |
| Microbial Cleanliness (Bioburden) | ≤30 cfu/g | | |
| Splash Resistance | Not Required | | Pass at ≥16.0kPa (≥120 mmHg) |

Summary:

| Performance Test as per EN 14683:2019 + AC:2019 | Savewo 3DMask Ultra |
|---|---------------------|
| | Type IIR |
| Bacterial Filtration Efficiency (BFE) % – <i>Staphylococcus aureus</i> (ATCC 6538) | Pass |
| Breathability (Differential pressure) | Pass |
| Microbial Cleanliness (Bioburden) | Pass |
| Splash Resistance | Pass |

Date : 2021-09-30
No. : HC21060355

Page 3 of 7

Test Result(s):

1. Bacterial Filtration Efficiency (BFE) %

Test method: EN 14683:2019 + AC:2019, Annex B

Summary: The BFE test is performed to determine the filtration efficiency of test articles by comparing the bacterial control counts upstream of the test article to the bacterial counts downstream. A suspension of *Staphylococcus aureus* was aerosolized using a nebulizer and delivered to the test article at a constant flow rate and fixed air pressure. The challenge delivery was maintained at $1.7 - 3.0 \times 10^3$ colony forming units (CFU) with a mean particle size (MPS) of $3.0 \pm 0.3\mu\text{m}$. The aerosols were drawn through a six-stage, viable particle, Andersen sampler for collection. This test method complies with EN 14683:2019 + AC:2019, Annex B

All test method acceptance criteria were met.

| Specimen(s) | Savewo 3DMask Ultra |
|-------------|---------------------|
| 1 | 99.6% |
| 2 | 99.6% |
| 3 | 99.7% |
| 4 | 99.8% |
| 5 | 99.9% |

- Notes :
- Challenge bacteria : *Staphylococcus aureus* (ATCC 6538)
 - Positive control average : 17.4 CFU
 - Negative control average : <1 CFU
 - Mean particle size : 3.1 μm
 - Testing side : Outside of specimen
 - Testing area : 49 cm²
 - Precondition : Minimum of 4 hours at (21 \pm 5) °C and (85 \pm 5) % relative humidity (RH)

Date : 2021-09-30
No. : HC21060355

Page 4 of 7

2. Differential Pressure

Test method: EN 14683:2019 + AC:2019, Annex C

Summary: The Differential Pressure test is performed to determine the breathability of test articles by measuring the differential air pressure on either side of the test article using a manometer, at a constant flow rate. This test complies with EN 14683:2019 + AC:2019, Annex C.

All test method acceptance criteria were met.

Sample : Savewo 3DMask Ultra

| Specimen(s) | Test area (in Pa/cm ²) | | | | | Average | |
|-------------|------------------------------------|------|------|------|------|--------------------|------------------------------------|
| | 1 | 2 | 3 | 4 | 5 | Pa/cm ² | mmH ₂ O/cm ² |
| 1 | 45.5 | 48.0 | 48.8 | 41.7 | 49.4 | 46.7 | 4.8 |
| 2 | 49.1 | 43.9 | 48.8 | 43.8 | 46.1 | 46.3 | 4.7 |
| 3 | 47.7 | 48.0 | 49.2 | 48.9 | 47.7 | 48.3 | 4.9 |
| 4 | 47.2 | 45.3 | 55.1 | 50.0 | 50.5 | 49.6 | 5.1 |
| 5 | 45.1 | 46.1 | 45.4 | 49.8 | 50.2 | 47.3 | 4.8 |

Notes : - 1 mmH₂O/cm² = 9.8 Pa/cm²
 - Flow rate: 8 Litre/min
 - Precondition : Minimum of 4 hours at (21±5) °C and (85±5) % relative humidity (RH)

Date : 2021-09-30
No. : HC21060355

Page 5 of 7

3. Microbial cleanliness (Bioburden)

Test method: EN 14683:2019 + AC:2019, Annex D and EN ISO 11737-1: 2018

Summary: The full mask is aseptically removed from the packaging and placed in a sterile 500ml bottle containing 300ml extraction liquid and shake for 5 min at 250 rpm. 100 ml of the extraction liquid is filtered through a 0.45µm filter and laid down on a TSA plate for the total viable aerobic microbial count. Another 100ml aliquot of the same extraction liquid is filtered in the same way and the filter plated on Sabouraud Dextrose agar (SDA) with chloramphenicol for fungi enumeration. The plates are incubated for 3 days at 30°C and 7 days at (20 to 25)°C for TSA and SDA plates respectively. The total bioburden is expressed by addition of the TSA and SDA counts. This test complies with EN14683:2019 + AC:2019, Annex D and EN ISO 11737-1: 2018.

All test method acceptance criteria were met.

Sample : Savewo 3DMask Ultra

| Specimen(s) | Mask Weight (g) | Total Bioburden (CFU/mask) | Total Bioburden (CFU/g) |
|--------------|-----------------|----------------------------|-------------------------|
| 1 | 4.13 | 0 | 0 |
| 2 | 4.20 | 0 | 0 |
| 3 | 4.13 | 0 | 0 |
| 4 | 4.16 | 0 | 0 |
| 5 | 4.17 | 0 | 0 |
| Mean: | | | 0 |

Notes : - CFU denotes Colony Forming Unit
- EN 14683:2019 + AC:2019 requirement: ≤ 30 CFU/g

Date : 2021-09-30
No. : HC21060355

Page 6 of 7

4. Synthetic Blood Penetration

Test method: ISO 22609:2004

Summary: This procedure was performed to evaluate surgical facemasks and other types of protective clothing materials designed to protect against fluid penetration. The purpose of this procedure is to simulate an arterial spray and evaluate the effectiveness of the test article in protecting the user from possible exposure to blood and other body fluids. The distance from the target area surface to the tip of the cannula is 30.5cm. A test volume of 2 mL of synthetic blood was employed using the targeting plate method.

This test method was designed to comply with ISO 22609:2004 (as referenced in EN14683:2019 + AC:2019) with the following exception: ISO 22609:2004 requires testing to be performed in an environment with a temperature of (21±5)°C and (85±5)% relative humidity (RH). Instead, testing was performed at ambient conditions within one minute of removal from the environmental chamber held at those parameters.

All test method acceptance criteria were met.

Test Pressure: ≥120mmHg

| Specimen Number | Sawewo 3DMask Ultra |
|--|---------------------|
| 1-20, 22, 24-26, 28-32 | None Seen |
| 21, 23, 27 | Seen |
| <p>Requirement: An acceptable quality limit of 4.0% is met for a normal single sampling plan when ≥ 29 of 32 test specimens show passing result (none seen)</p> | |

Notes : - Test Side: Outside
- Precondition : Minimum of 4 hours at (21±5) °C and (85±5) % relative humidity (RH)


STC

Test Report

Date : 2021-09-30
No. : HC21060355

Page 7 of 7

Photo(s):

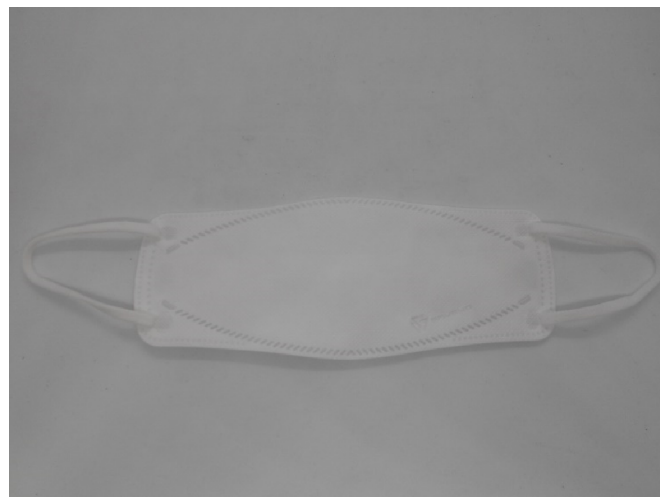

***** End of Test Report *****

Conditions of Issuance of Test Reports

1. All samples and goods are accepted by The Hong Kong Standards & Testing Centre Limited (the "Company") solely for testing and reporting in accordance with the following terms and conditions. The Company provides its services on the basis that such terms and conditions constitute express agreement between the Company and any person, firm or company requesting its services (the "Clients").
2. Any report issued by the Company as a result of this application for testing service (the "Report") shall be issued in confidence to the Clients and the Report will be strictly treated as such by the Company. It may not be reproduced either in its entirety or in part and it may not be used for advertising or other unauthorized purposes without the written consent of the Company. The Clients to whom the Report is issued may, however, show or send it, or a certified copy thereof prepared by the Company to his customer, supplier or other persons directly concerned. Subject to clause 3, the Company will not, without the consent of the Clients, enter into any discussion or correspondence with any third party concerning the contents of the Report, unless required by the relevant governmental authorities, laws or court orders.
3. The Company shall be at liberty to disclose the testing-related documents and/or files anytime to any third-party accreditation and/or recognition bodies for audit or other related purposes. No liabilities whatsoever shall attach to the Company's act of disclosure.
4. The Company shall not be called or be liable to be called to give evidence or testimony on the Report in a court of law without its prior written consent, unless required by the relevant governmental authorities, laws or court orders.
5. The results in Report apply only to the sample as received and do not apply to the bulk, unless the sampling has been carried out by the Company and is stated as such in the Report.
6. When a statement of conformity to a specification or standard is provided, the ILAC-G8 Guidance document (and/or IEC Guide 115 in the electrotechnical sector) will be adopted as a decision rule for the determination of conformity unless it is inherent in the requested specification or standard, or otherwise specified in the Report.
7. In the event of the improper use the report as determined by the Company, the Company reserves the right to withdraw it, and to adopt any other additional remedies which may be appropriate.
8. Sample submitted for testing are accepted on the understanding that the Report issued cannot form the basis of, or be the instrument for, any legal action against the Company.
9. The Company will not be liable for or accept responsibility for any loss or damage howsoever arising from the use of information contained in any of its Reports or in any communication whatsoever about its said tests or investigations.
10. Clients wishing to use the Report in court proceedings or arbitration shall inform the Company to that effect prior to submitting the sample for testing.
11. Subject to the variable length of retention time for test data and report stored hereinto as to otherwise specifically required by individual accreditation authorities, the Company will only keep the supporting test data and information of this test report for a period of three years. The data and information will be disposed of after the aforementioned retention period has elapsed. Under no circumstances shall we provide any data and information which has been disposed of after the retention period. Under no circumstances shall we be liable for damages of any kind, including (but not limited to) compensatory damages, lost profits, lost data, or any form of special, incidental, indirect, consequential or punitive damages of any kind, whether based on breach of contract of warranty, tort (including negligence), product liability or otherwise, even if we are informed in advance of the possibility of such damages.
12. Issuance records of the Report are available on the internet at www.stc.group. Further enquiry of validity or verification of the Reports should be addressed to the Company.