



Test Report

No.T32120280386SN

Date: Aug 02, 2021

Page 1 of 3

SAVEWO LIMITED

1/F., 266-270 TEXACO ROAD, TSUEN WAN, N.T., HONG KONG

The following samples were submitted and identified by/on behalf of the client as:

3D MEOW KIDULTS, 3D MEOW S, 3D MEOW L, 3D BEAR

Case No. : CA321202821099
Lot No. / Batch Code : NOT PROVIDED
Sample Description : WHITE FABRIC
Style / Item No. : 3D MEOW-KDT, 3D MEOW-S2, 3D MEOW-L2, 3D BEAR BB
Manufacturer : SAVEWO LIMITED
Country of Origin : HONG KONG
Sample Receiving Date : JUL 15, 2021
Testing Period : JUL 15, 2021 – JUL 30, 2021

Test Requested	Conclusion
Viral filtration efficiency (VFE)	See Result

***** FOR FURTHER DETAILS, PLEASE REFER TO THE FOLLOWING PAGE(S) *****

Signed for and on behalf of
SGS Hong Kong Ltd.

Au Kam Chi, Gigi
Technical Manager

Signed for and on behalf of
SGS Hong Kong Ltd.

Tsang Chuk Hai
Senior Microbiologist

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.



Viral filtration efficiency (VFE)

Test Side : Inside
 Test Method : With reference to modified method of
 ASTM F2101-19 / EN 14683:2019+AC:2019 Annex B
 Pre-Conditioning : Minimum of 4 hours at 21±5°C and 85±5% R.H.
 Dimensions of test specimen : 100 mm x 100 mm
 VFE Test Area : 49 cm²
 VFE Flow Rate : 28.3 l/min
 Test Virus : Escherichia coli bacteriophage Phi X174 (ATCC 13706-B1)
 Mean Particle Size : 2.8 µm
 Positive Control Average : 1.1 x 10³ CFU
 Negative Monitor Count : < 1 CFU

Test Specimen	Percent VFE (%)
1	99.9
2	99.9
3	99.9
4	99.9
5	99.9


Note : Please note count total for each stage can be provided upon request

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Sample Photo:

Sample Picture (As received)

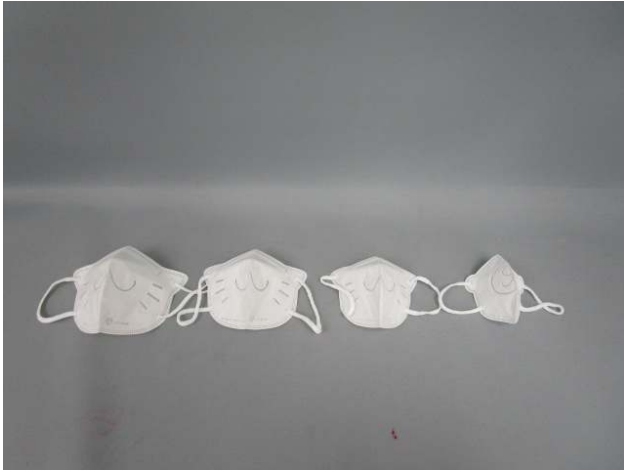

SGS authenticate the photo on original report only

*** End of Report ***

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.