

Test Report

No.T32120280391SN

Date: Aug 11, 2021

Page 1 of 3

SAVEWO LIMITED

1/F., 266-270 TEXACO ROAD, TSUEN WAN, N.T., HONG KONG

The following samples were submitted and identified by/on behalf of the client as:

3D MEOW L

Case No. : CA321202821069
 Lot No. / Batch Code : NOT PROVIDED
 Sample Description : WHITE MASK
 Style / Item No. : 3D MEOW-L2
 Manufacturer : SAVEWO LIMITED
 Country of Origin : HONG KONG
 Sample Receiving Date : JUL 15, 2021
 Testing Period : JUL 15, 2021 – AUG 03, 2021

Test Requested	Conclusion
Microbial cleanliness (Bioburden)(EN 14683:2019+AC:2019 Annex D)	PASS

***** FOR FURTHER DETAILS, PLEASE REFER TO THE FOLLOWING PAGE(S) *****


Signed for and on behalf of
 SGS Hong Kong Ltd.



Wong Kin Man, Gilman
 Technical Development Manager

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Document.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

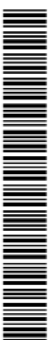

Test Report

No.T32120280391SN

Date: Aug 11, 2021

Page 2 of 3

Test Results:

Microbial cleanliness (Bioburden)(EN 14683:2019+AC:2019 Annex D)

Test Methods

Bioburden

The analyses were performed according to EN 14683:2019+AC:2019 and ISO 11737-1:2018

Test Results

SGS Sample No.:HKHC210700002188-101

Article Number	Mask Weight	Total Bioburden, cfu/mask	Total Bioburden, cfu/g
1	2.00g	12	6.00
2	2.10g	18	8.57
3	2.07g	18	8.70
4	2.07g	42	20.29
5	2.07g	21	10.14
Mean:		22.2	10.7

Recovery Efficiency	Correction Factor
45.2%	2.2

Microbial Cleanliness (Bioburden): 23.7 cfu/g

Standard requirement#: ≤ 30 cfu/g

Note:


1. Results reported on the submitted sample on an as received basis.
2. cfu = Colony Forming Units
3. Extraction method: by shaking at 250rpm for 5 minutes
4. # EN 14683:2019+AC:2019 - Medical face masks - Requirements and test methods – Performance requirements for medical face masks – Microbial cleanliness

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Document.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Sample Photo:

Sample Picture (As received)

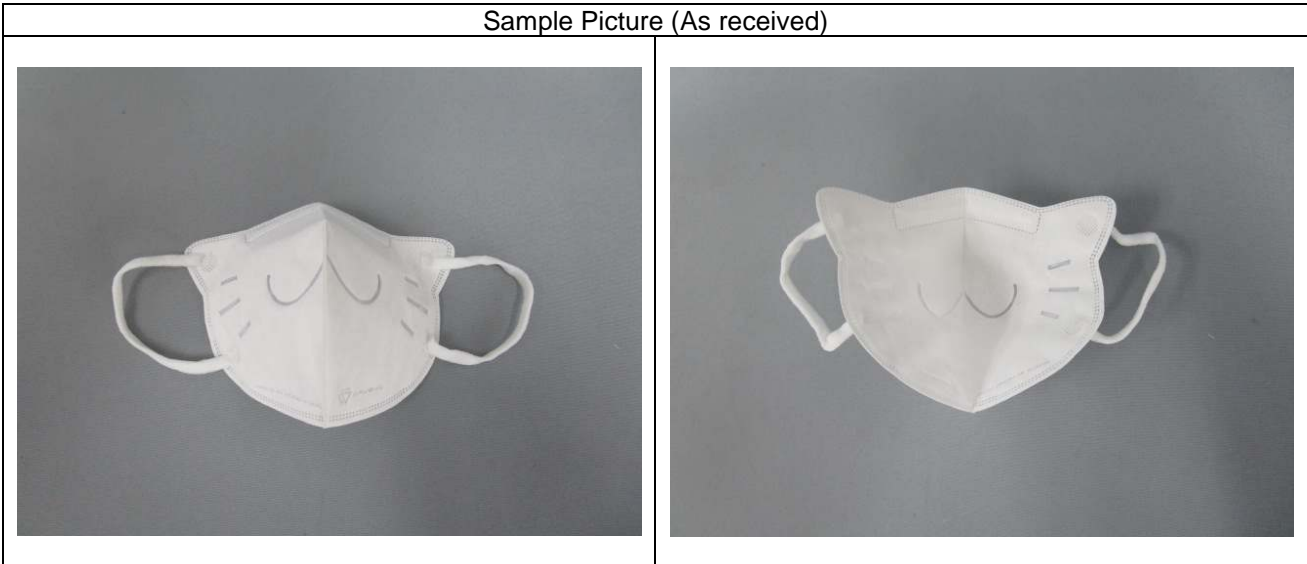

SGS authenticate the photo on original report only

*** End of Report ***

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Document.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.